


OXYGEN PLANT

O₂-JEN



How Is Work?

O2-jen plant, produce high purity oxygen on 5 bar and supply your gas system.

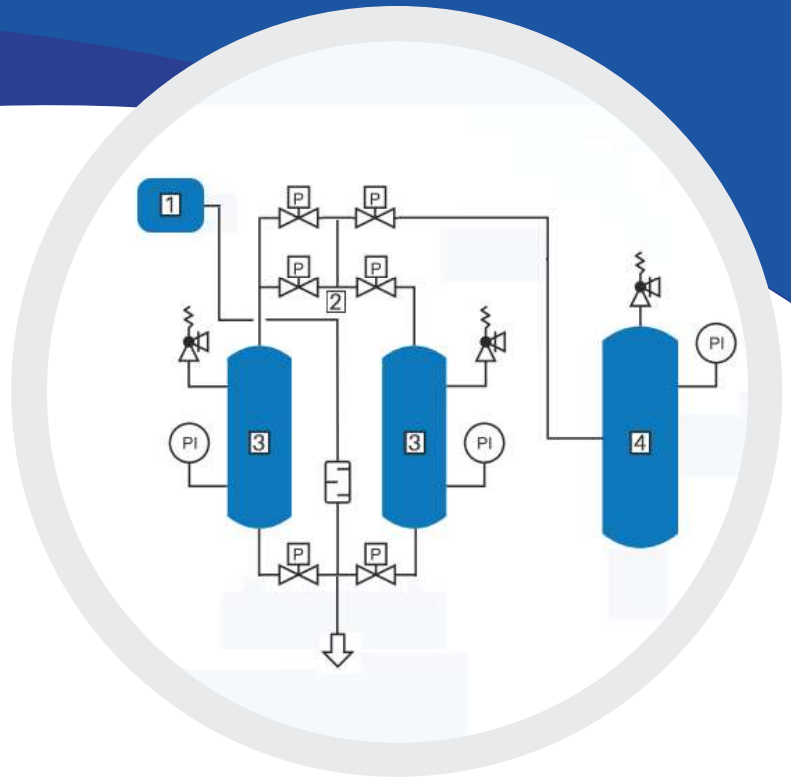

Air compressor (1) supply the **O2-jen plant** with air , **O2-jen plant (3)** has got zeolite inside with this, keep O₂ and separate nitrogen. The valve group (2) send the O₂ to oxygen tank (4) and exit the nitrogen outside. In oxygen tank always 5 bar oxygen on %95 purity ready to use.


WHY DO YOU PAY MORE?

When you are use oxygen cylinder or liquid oxygen , you dont pay these gas, on same time

- On cylinders there are high pressure so you get this high pressure risk
- When you change cylinder there can be change pressure on pipeline
- You rent cylinders
- You pay delivering charges
- You lost time and people power when change cylinder
- You pay every months on oxygen, every time you lost time for tenders



FORGET CYLINDER TENDERS AND COSTS

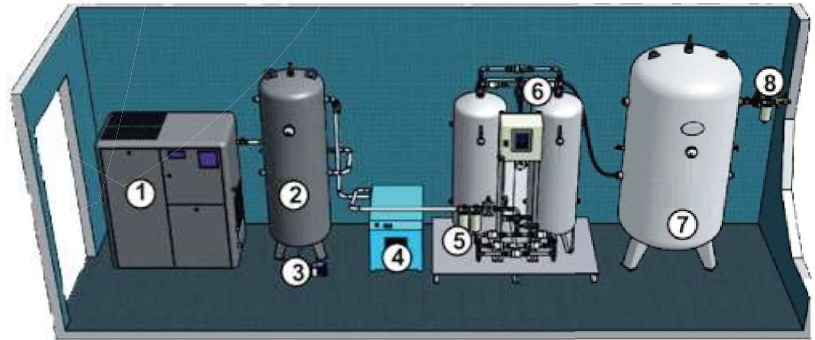

Medikal Oxygen produce

SCREW TYPE AIR COMPRESSOR

1

AIR TANK

2

AUTOMATIC WATER DISCHARGE

3

AIR DRYER

4

AIR FILTER

5

O2-jen plant

6

O2 tank

7

Exit Filter

8

- Working in PSA
- Automatic and manual working modes
- Easy installation and use
- %95 pruity



Performance

Purity(%)	%94
Oxygen flow (nm3/hr)	4-50

PRODUCE YOUR OWN OXYGEN

ABOUT US

O₂-JEN

 **RAVZA**
MEDİKAL

*YOU CAN SAVE YOUR
TIME AND MONEY WITH
O₂-JEN PLANT . YOU
CAN PRODUCE YOUR
OWN OXYGEN AND YOU
WILL FORGET OXYGEN
CYCLINDERS*

 +90 506 688 40 03

 murat@ravzamedikal.com.tr

**You can contact us
for free pattern de-
tection**